


1. Making a template from an existing page

- Press the 'Tools' icon button
- Select 'create page'
- Give the page a name e.g. Template 36
- Press OK
- Select 'copy from an existing page' (you can only copy a page from your own User Area)
- Select the 36 button page you wish to copy as your template
- Modify the page to customise your template e.g. keep the function buttons like a home 'page link' button and a clear 'tool' button, but clear the vocabulary.
- This template is now ready to use whenever you make a new page (create page, name it, copy from an existing page, select your template)

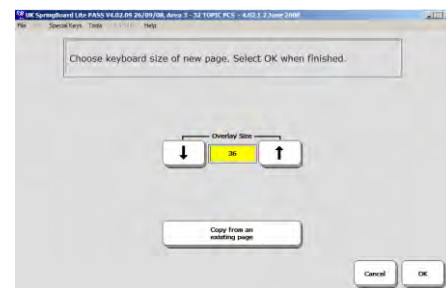

2. Making a template from scratch

- Press the 'Tools' icon button
- Select 'create page'
- Give the page a name e.g. Template 8
- Press OK
- Select overlay size e.g. 8 buttons
- Select OK
- Modify your page, by adding e.g. a home 'page link' button (so you can get out of the page), clear 'tool' button, change colour scheme etc.
- You now have a customised 'blank template' to use. You could also make templates of 4, 15 or 36 buttons to use within your User Area (copy from an existing page, select the template)